

## PIKA FASTORANGE® B TUBES 48-PACK

Incubation tubes for the detection of beer spoiling microorganisms

Cat. No. 2036-10

Description	Amount	Storage
Incubation tubes for the detection of beer spoiling bacteria and Dekkera (Brettanomyces) yeast.	48 x 5 mL tubes	Store dark at room temperature

**Warning!** Read the manual and the Safety Data Sheets before starting the analysis. Safety Data Sheets are available in the download area from [www.pika-weihenstephan.com](http://www.pika-weihenstephan.com). All handling steps should be performed under sterile conditions. Wear appropriate protective clothing  
For *in vitro* use only.


### Product description

PIKA FastOrange® B Tubes are incubation tubes prefilled with culture medium. They are easy to use and directly applicable at the sampling point.

PIKA FastOrange® B Bouillon is a culture medium developed to detect microbial contamination during the brewing process and in brewery products.

Acid producing beer spoiling bacteria as lactic acid bacteria and Pediococci are easily detected by a color change of the culture medium from violet to yellow.

In case no other viable microorganisms are present in the sample, also Dekkera (Brettanomyces) yeasts may be detected. For the specific detection of Brettanomyces yeasts we recommend FastOrange® BRETT Tubes (Cat. No. 2037-10).

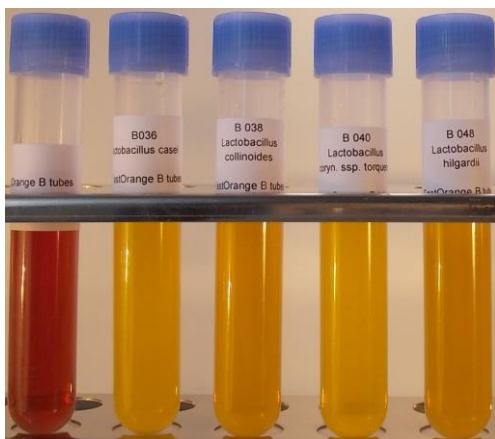

### Detectable microorganisms

Microorganism	Growth conditions
Lactic acid bacteria (Lactobacillaceae) - Lactobacillus sp. - Pediococcus sp.	aerobic and anaerobic
Acetic acid bacteria (Acetobacteriaceae) - Acetobacter sp. - Gluconobacter sp.	aerobic
Gram negative beer spoiling bacteria, including - Pectinatus sp. - Megasphaera sp.	anaerobic
Beer spoiling yeast, including - Dekkera sp. (Brettanomyces)	aerobic and anaerobic

Growth of brewer's yeast and most other yeasts is suppressed.

### Guidelines for use

- Under sterile conditions, pour liquid directly into the tube with medium and mix.  
- 3-8 mL liquid sample volume is recommended
- Incubate tubes at  $25 \pm 2^\circ\text{C}$  for 2-7 days.

**Important!** For the detection of anaerobic microorganisms, follow one of the following recommendations

- Fill the tube up to the brim and then close it
- Alternatively, fix the cap slightly on the tube to allow gas exchange and then incubate tube in a suitable anaerobic vessel.

### Results of visual evaluation

Target microorganism	Samples have to be regarded as positive if:
Acid producing bacteria including beer spoiling bacteria	color change: violet-brown to yellow
Other microorganisms including indicator- and accompanying flora	turbidity and potential sediment formation

Growth of beer spoiling microorganisms is indicated by turbidity together with color change, and often additionally by sediment formation.

### We recommend

- Microscopic examination and/or PCR analysis to verify the presence of beer spoiling bacteria in positive enrichments.
- Serial enrichment if an immediate color change occurs after the sample is mixed with broth.

Serial enrichment	Method
Second enrichment	1. After 3 days of enrichment, transfer half the volume of the enriched sample into a fresh Fast Orange® B Tube 2. Incubate at $25 \pm 2^\circ\text{C}$ for 3-5 more days

## General information

Store the product in the dark at room temperature (max. 25°C). Cooling below 25°C is NOT necessary.

Due to manufacturing, slightly differences in color of culture medium may occur within lots. This is NOT influencing product quality.

Best before date for unopened product is given on the outer label. After opening we cannot guarantee for shelf life.

The product is not suitable for human or animal consumption. It must not be used for the direct propagation of microorganisms which later are used for food production or might get into contact with food.

## FastOrange® B Products

B Bouillon	(12 x 240 mL)	Item #2036-1
B Agar	(12 x 170 mL)	Item #2036-2
B Ready to Use Tubes 48-pack	(48 x 5 mL) with 48 Swabs	Item #2036-3
B Ready to Use Tubes 24-pack	(24 x 5 mL) with 24 Swabs	Item #2036-14
B Tubes 48-pack	(48 x 5 mL)	Item #2036-10
B Tubes 24-pack	(24 x 5 mL)	Item #2036-15
B Enrichment Bottles	(15 x 50 mL)	Item #2036-11



PIKA Weihenstephan GmbH  
Raiffeisenstrasse 31A  
85276 Pfaffenhofen  
GERMANY  
Phone +49 (0) 8441 879 48 30  
Fax +49 (0) 8441 879 48 31  
[www.pika-weihenstephan.com](http://www.pika-weihenstephan.com)  
[order@pika-weihenstephan.de](mailto:order@pika-weihenstephan.de)

---

**Notes:** The relevant antibiotics/fungicides contained in the medium fall short of critical values that require monitoring or declaration according to regulation (EG) 1907/2006 (REACH). When properly applied, the medium may be disposed of through the normal sewage system. It is strongly recommended to inactivate the live microorganisms in any enriched sample by heating to 121°C/250°F for 20 min (autoclave) to exclude a release of live microorganisms. Although this information was collected thoroughly, we cannot guarantee that any of the content is incomplete or incorrect. We do not take over any warranty for consequences which are resulting from improper handling or incorrect use of this product. Additionally, always comply with the applicable laws, regulations and directives of your country. PIKA Weihenstephan® and FastOrange® are registered trademarks in Germany and other countries.